

Power Tower View – Complete practice supply
with a minimal footprint.



Power Tower View – The perfect solution for every requirement

Compressed air and suction combined with a minimal footprint.

The compact and quiet Power Tower View can be installed flexibly in the surgery with a footprint of just 0.4 m², saving valuable space over the installation of equivalent separate units – a real advantage given the expense of current rental costs. This unit is also the perfect fit for any surgery with its high-quality design and low noise level thanks to the sound insulation. This leaves you free to position the unit anywhere you want in the surgery. Thanks to the many individual configuration options, Power Tower View is a flexible solution that will grow with the surgery and adapt effortlessly to new circumstances.



Compressed air at the top, suction at the bottom

At the top is where the heart of the surgery beats: this is where the powerful and durable compressor reliably delivers oil-free, dry and hygienic compressed air. The air quality of our compressors in the Power Tower View complies with the requirements of the new ISO 22052 for dental compressed air. The effective Dürr Dental suction systems – in either dry or wet variants – for continuous and ergonomic treatments safe from the risk of infection are housed at the bottom. The Power Tower View is also certified as a Class IIa medical device in accordance with the MDR (EU 2017/745).



This is what counts:

- Complete practice supply with a minimal footprint thanks to the compact, modern and appealing design that will help your surgery rooms look even better
- Increased operational reliability thanks to early warning messages and automatic maintenance reminders (perfectly maintained units, fast troubleshooting, remote servicing)
- Demand-orientated compressor control offers redundant operation
- User-friendly adjustment of device settings to the conditions in the practice
- Soundproofing enclosure keeps noise to a comfortable level of max. 56 dB(A)
- Variable configuration options make the unit perfect for any surgery
- User-friendly and comfortable operation of the unit



A solution with options for more

The modular concept of Power Tower View offers flexible options for upgrading and modification. For example, you can retrofit an additional compressor unit in order to boost output or create redundancy for increased operational reliability.

One for all – and everything in a single unit

Power Tower View combines a compressor unit, the compressed air tank, the suction unit, the membrane drying unit and – if required – also the amalgam separator into a single, compact unit offering perfect thermal and physical conditions. The modular concept will fulfil numerous individual requirements and allows the system to be set up for up to 4 therapists, providing the perfect solution for any size of surgery.

Stay relaxed with a clear view of everything – with VistaSoft Monitor



Our software integrates the supply systems in the surgery network, ensuring that a constant overview of the current status of monitored systems is always provided. The surgery is also automatically notified about upcoming maintenance tasks, such as filter changes on the compressor or when the amalgam container reaches a fill level of 95%. In the event of a fault, the software will indicate whether you can remedy the problem yourself, or it offers the option of informing a service technician directly via e-mail. The technician can immediately look for a solution using a remote maintenance function.



Scan me



Technical data at a glance

	VS 600	VS 900 S	VS 900 S	VS 1200 S	V 600	V 900 S	V 1200 S
Number of users	2	3	3	4	2	3	4
Compressor	Duo	Duo	Duo	Duo	Duo	Duo	Duo
2nd unit (optional)	✓	✓	✓	✓	✓	✓	✓
Amalgam separator (optional)	✓	✓	✓	✓	X	X	X
Delivery at 5 bar [l/min] with 1 unit	115	115	115	115	115	115	115
Delivery at 5 bar [l/min] with 2 units	230	230	230	230	230	230	230
Noise level [dB(A)] at an ambient temperature of 35°C	56	56	56	56	56	56	56
Voltage (V)	230	230	400	400	230	400	400
Frequency (Hz)	50	50	50	50	50	50	50
Dimensions (H x W x D cm)	208 x 64 x 61	208 x 64 x 61	208 x 64 x 61	208 x 64 x 61	208 x 64 x 61	208 x 64 x 61	208 x 64 x 61
Weight (kg)	275	275	275	275	265	265	265

* Delivery measured without dry air



GERMAN
DESIGN
AWARD
SPECIAL



GERMAN
INNOVATION
AWARD '21
WINNER



Accessories


Filter finder overview website: <https://www.duerredental.com/en/products/compressed-air/materials-and-accessories/filter-finder/>



Air intake filter for
compressor units



Bacteria filter for
membrane drying unit

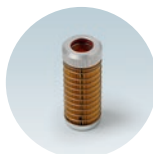

Fine filter for
membrane drying unit



Coalescence filter for
membrane drying unit

DÜRR DENTAL SE
Höpfigheimer Str. 17
74321 Bietigheim-Bissingen
Germany
More information under www.duerredental.com
info@duerrdental.com



Scan me

**DÜRR
DENTAL**
THE BEST, BY DESIGN